

### 6 IN 1 HIFU & LIPOHIFU SYSTEM

HIFU


Face Lift   Wrinkle F	Removal   Body S	Slimming
-----------------------	------------------	----------



### How lipohifu works?

#### Introduction

The Lipo hifu system uses high-intensity focused ultrasound deliver uniform energy at precise depth (13mm) to destroy targeted subcutaneous fat non-invasively, causing a rapid local temperature rise, destroying targeted adipose tissue just beneath the skin in the abdomen and flanks, predictable results after a single one-hour treatment. It can penetrate the skin layers and reach the targeted fat tissue without harming skin or surrounding tissues. Meanwhile the metabolized contents (triglyceride, fatty acid) of cells will excrete out of body automatically by blood circulation and lymphatic drainage without harming human body.

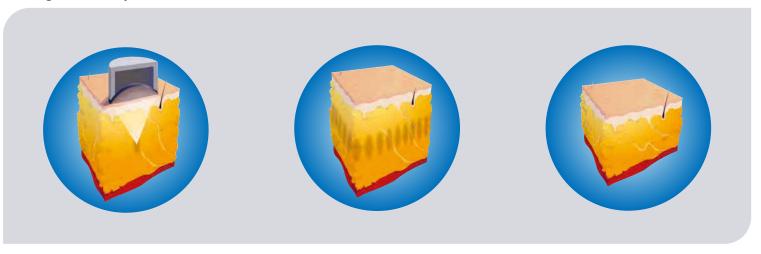

### **How HIFU Skin Tightening works?**

#### Introduction

High Intensity Focused Ultrasound (HIFU) directly delivers heat energy to skin and subcutaneous tissue that can stimulate and renew the skin's collagen and thus consequently improving the texture and reducing sagging of the skin. It literally achieves the results of a facelift or a body lifts without any invasive surgery or injections, moreover, an added bonus of this procedure is that there is no downtime.


HIFU(High Intensity Focused Ultrasound) is delivered into the deep skin fractionally



Cause thermal coagulation to SMAS (Superficial Muscular Aponeurotic System), fascia, deep dermis and subcutaneous fibrous tissue.



By wound healing, skin tightening and rejuvenation happens gradually with time

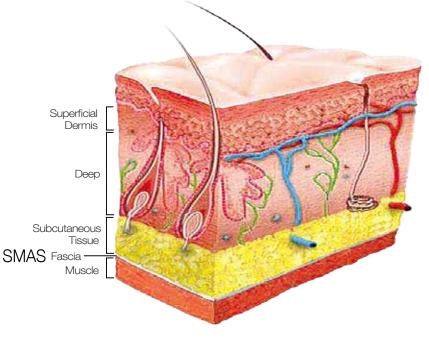

### 12D HIFU

20,000 SHOTS

**HIFU** cartridges for

body and face










3.0mm 4.5mm 9.0mm 12.0mm 6.0mm

Single line

Multi line

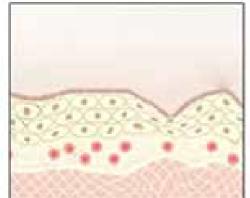

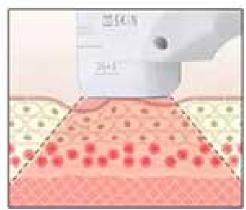

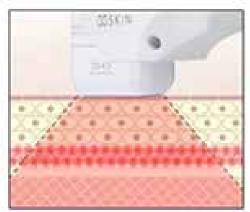



### 600 SHOTS

**Lipo HIFU cartridges** 

fat loss and body slimming










**LIPOHIFU** 

**HANDLE** 

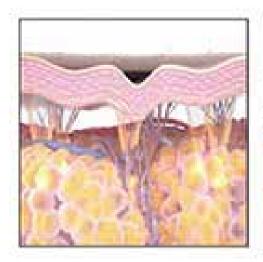

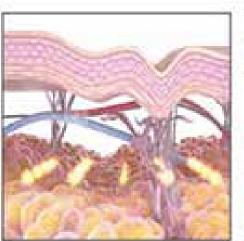

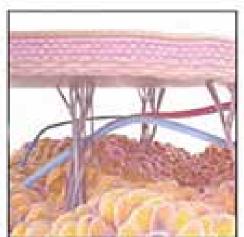

# 360° FACE HIFU CARTRIDGE

### 40,000 SHOTS

**HIFU** cartridges for

Face lift&compact



1.5mm

3.0mm

## DEP INJECTION HANDLE

Vibration mode to increase blood circulation and comfort during operation



Non injection, noninvasive treatment



Effective introduction of macromolecular components



Instant introduction, no harm to the skin



Vibration mode to increase blood circulation and comfort during operation




### FRACTIONAL RF HANDLE LED Orientation



Different tips targeted to different body part.


To create the most copact effect.

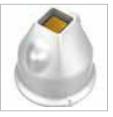

S for eyes (5\*5mm)



M for face (15x15mm)

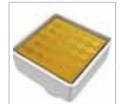

L for body (34.5x34.5mm)




The patented Flow Needling Technology serially inserts the microneedles row-by-row for more comfortable procedure.



2. Serial insertion reduces unnatural stretching in order to further reduce pain and downtimes. The depth of the coagulation zone is precisely controlled.




 Sub-dermal delivery of RF energy thermal damage to the epidermis providing faster healing time for your patients.

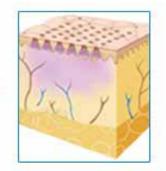

4. The adjustable depth of the RF energy delivery allows for multiple passes and treatment of delicate areas.

### VAGINAL EXTENSION HANDLE

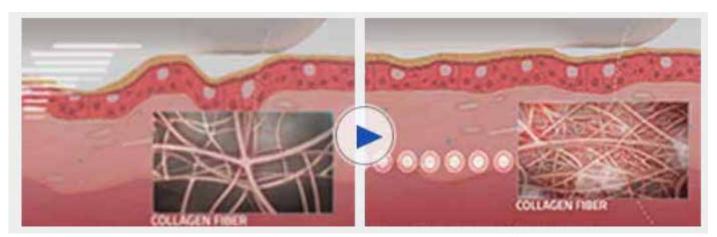

Auto-decontrolled rotation angle of 360









before

In operation

after






### Treatment range

Is a safe an effective means of non-invasive body sculpting.

- ✓ Parts all over the body
- Dissolve abdominal fat, reduce abdominal fat prominent
- Face, neck, arms and legs, buttocks and abdomen
- Postpartum and after liposuction reshape.
- Slimming, firming, improve cellulite and stretch marks

### FACIAL AND BODY TREATMENT



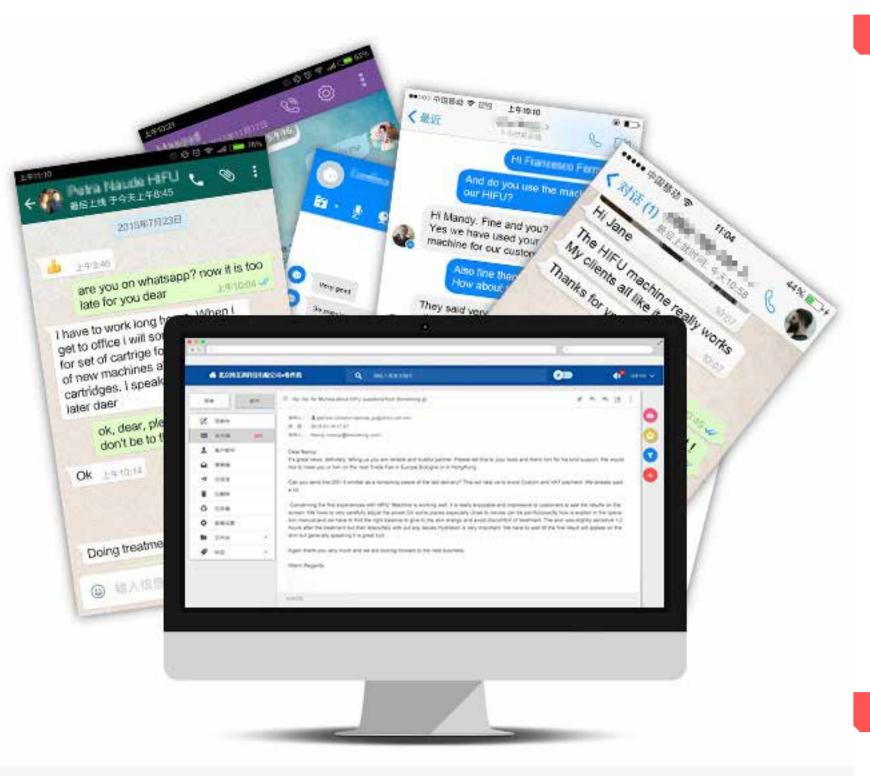

Concerning the first experiences with FU5. Machine is working well. it is really enjoyable and impressive to customers to see the results on the screen.

### Body parts ▼



#### Face ▼



### HIFU & LIPOHIFU SYSTEM Non-Invasive Fat Reduction



